


# For All Breeders, Kennels and Vets

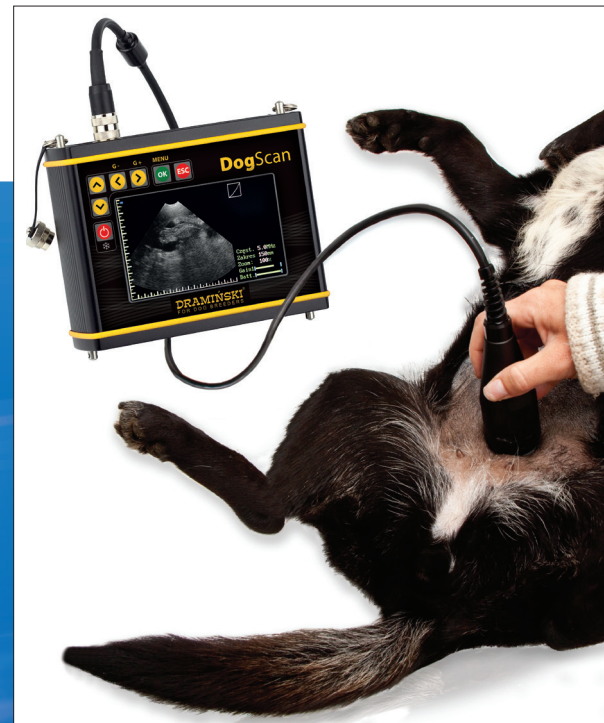

## DogScan

A portable LCD Ultrasound scanner to quickly and accurately diagnose pregnancy in small animals.

### Application:


- This small, lightweight and robust scanner has been specifically designed for dog breeders to enable them to examine their own bitches, to confirm or exclude pregnancy.

### Benefits:


- Early pregnancy detection,
- Monitoring pregnancy development,
- Professional care of pregnant bitch.

### Equipment Features:

- Small, light and handy – 1300 g only,
- Robust and durable construction,
- Shock-proof and humidity resistant aluminium casing,
- Touch sensitive membrane keyboard,
- 5.0", LCD LED screen showing high quality clear images,
- Detached mechanical sector 5MHz abdominal probe,
- External rechargeable battery with 6 hours of continual usage,
- Charging time 2,5 h.

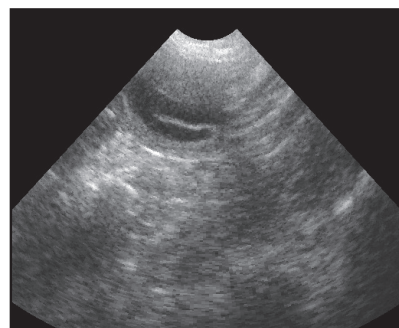

Day 25th of pregnancy

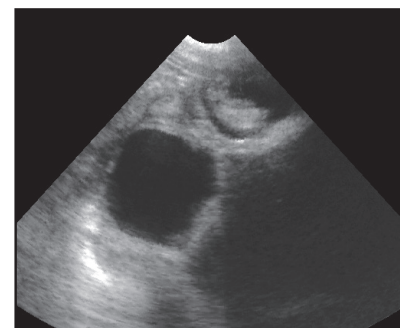

Day 28th of pregnancy



Day 36th of pregnancy







## DOD2 Dog Ovulation Detector

Helps to determine the best mating time, reduces time and travel costs, research and repeated mating, identifies the so-called „silent oestrus“ and reproductive process disorders, eliminates the stress associated with repeated mating attempts.

### USB interface

- remote software update and data transfer to a computer

### Two operating modes

- standard / advanced.

### Easy to disinfect

- washed under running water, which makes it easier to keep clean.

### The kit includes:

- DOD2 device with lanyard,
- carrying case with a foam insert,
- manual,
- 4 AA-type batteries 1.5 V (LR6),
- USB cable for communication with a computer.

### Main characteristics:

- reliable results,
- large display with energy-saving LED backlighting,
- safe for animals and humans,
- does not require periodic calibration,
- memory for 100 dogs and 200,000 measurements with date and time,
- data transmission to a computer for easy interpretation of results (oestrus graphs, printouts, reports, archiving).
- 2-year warranty,

### Additional features:

- real-time clock,
- pop-up menu,
- setting the display operating parameters,
- language selection,
- automatic shutdown time change,
- measurement settings.



Free software available for download from manufacturers website.



## Pregnancy Detector for Dogs


A valuable management and diagnostic tool in modern breeding. This portable detector is a small, easy to operate device, powered by a common 9 volt PP3 battery, providing a very fast and reliable method of testing bitches for pregnancy. The waves emitted by the probe, which is placed on the skin in the position as shown on the diagram, are reflected by the uterine wall of the pregnant bitch. When testing, the signals are then returned to the unit, analysed and when pregnancy is detected, both the diode and acoustic signal rate increases rapidly.

### The Detector:

- Identifies bitches that are not pregnant,
- Helps to identify bitches with problems in reproduction.

### Benefits:

- Helps to confirm or exclude pregnancy in bitches,
- Confirms efficiency of insemination,
- Supports a breeder in decision making,
- Manages and reduces cost of animal care.

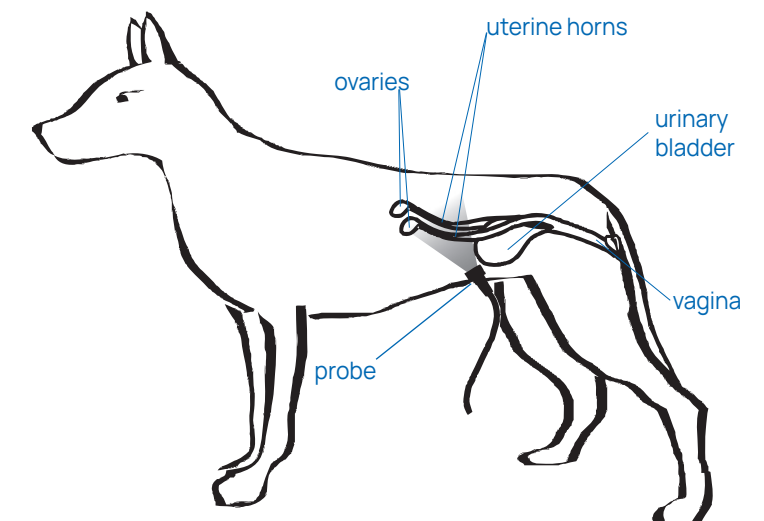

The figure shows the correct application of the probe to the bitch's stomach.



## 4Vet Slim

Ultrasound scanner with a high quality imaging for the diagnosis of small animals and horses. Perfect for veterinary mixed practises.


### Benefits:

- Obtain a clear image no matter where you are, thanks to a high contrast, bright screen.
- Perform examination even with limited access to electricity.
- All the necessary functions are at hand, a sensitive touch screen makes setting of a proper image faster,
- After 25 seconds Slim is ready for operation,
- Wide choice of the probes enables examination of different patients,
- Two-year warranty, for both the device and the probes.


### Equipment Features:

- Touch screen 10,4",
- Ready to work in 25 sec.,
- Compatible with various probes and OPU set,
- Needle enhancement,

- Color Doppler,
- Modern design,
- AC plug or battery powered,
- Up to 4 hours of continuous work on battery.

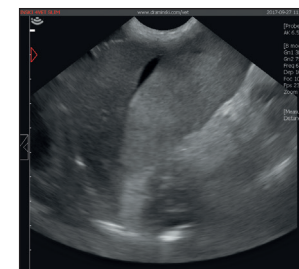

Free liquid in abdomen

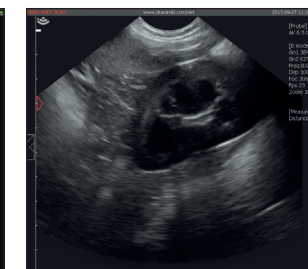

Kidney of healthy dog



Adrenal gland of healthy dog

