

iScan

ULTRASOUND SCANNER DESIGNED FOR VETS WORKING IN EXTREME FIELD CONDITIONS

- perfect image quality
- operates with a fixed electronic linear probe
- small, light and portable
- extremely efficient batteries
- ready to work immediately
- waterproof



DRAMIŃSKI *iScan*

ULTRASOUND SCANNER DESIGNED FOR VETS
WORKING IN EXTREME FIELD CONDITIONS



- ◀ It has perfect image quality guaranteeing precision and reliability in testing cows and mares.
- ◀ It operates with a fixed electronic linear probe.
- ◀ It is small, light and portable, which guarantees comfort of operation.
- ◀ It has extremely efficient batteries enabling 5 hours of continuous work.
- ◀ Once taken out of the case, it is ready to work immediately.
- ◀ It is waterproof, it can be easily washed under a running tap.

NOTE: the device is also available for left-handed people.



PORTABLE ULTRASOUND SCANNERS TO WORK IN THE FIELD

Standard package includes:

ultrasound scanner with a fixed linear probe, 1 external battery, battery charger, USB cable, carry strap, transport case, software to transmit images to PC, user manual.

Optional package: ultrasound gel (250 ml), sun protection shield, gel pad to test tendons and eyeballs at horses, additional external battery.

Technical specification:

- Size: 17.0 x 17.0 x 6.0 cm
- Weight: 1720 g with a probe
- 5" LED-backlit LCD display

Image:


- B, B+B, B+M,
- zoom 100 – 160%,
- freeze,
- regulation: gain, gamma, frequency, focus.

Probes:


- linear endorectal 4 – 9 MHz, penetration depth up to 12 cm

Functions:

- save images and cine loop with description and date to the hard disc,
- transmission of images and cine loop to PC,
- measuring (distance, surface, volume, fetal biometry charts),
- regulation of scanning range,
- presets.

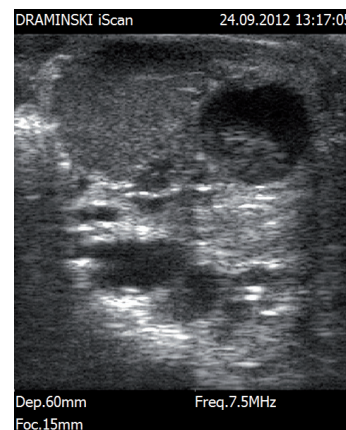

Cow – ovary

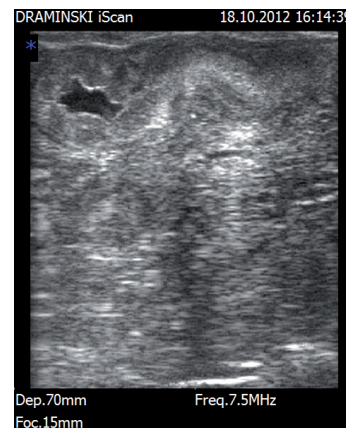

Heifer – the 28th day of pregnancy

DRAMIŃSKI S.A.

Wiktora Steffena 21, 11-036 Sząbruk, Poland
tel. +48 89 675 26 00, sales@draminski.com

**WE CONSTRUCT
PRODUCE
SELL
SERVICE**